

QUICKIE[®]

QM-7 Series

Full-featured mid-wheel drive performance
with proven reliability.



Live Without Limits™

The Quickie QM-7 series offers a high-performing, reliable, and stylish mid-wheel drive power solution for a client who qualifies for a Group 3 Standard (QM-710), Group 3 Heavy Duty (QM-715 HD) or Group 4 High Speed (QM-720) power mobility base. Features like SpiderTrac™ suspension, an intuitive mid-wheel drive base and PG Drives electronics provide superior ride comfort over any terrain. The QM-7 series is available with a full range of sophisticated no-power option, single-power option and multi-power option seating systems designed for complex clinical power seating applications. The result: a comprehensive modular approach that combines performance with versatility.



Power Base— Designed for Performance

- Mid-wheel drive intuitive handling
- Linix 4-pole standard motor package for all models
- Contoured, sleek base style available with 5 shroud color options



- Independent 6-wheel suspension for superior traction and handling across multiple terrains
- Gas dampening system for smooth, stable transitions
- 3" curb climbing capability
- 22° transition angle helps prevent high centering on curb cuts

Available in a wide range of models to meet requirements of many Group 3 and Group 4 codes:

GROUP 3: Standard Wt. Cap. (6 mph, 300 lbs.)

- Quickie QM-710 BC:
No power Option Base w/ ASAP II Rehab Seat
- Quickie QM-710 SC:
Single Power Option Tilt (SPOT) or Recline base with ASAP II or Recline Rehab Seat
- Quickie QM-710 MPC:
Multiple Power Option Base (SPOT Tilt and Recline) with Recline Rehab Seat

GROUP 3: Heavy Duty (5 mph, 400 lbs.)

- Quickie QM-715 HD BC:
No power Option Base w/ ASAP II Rehab Seat
- Quickie QM-715 HD SC:
Single Power Option Tilt (SPOT) with ASAP II

GROUP 4: High Activity (8 mph, 300 lbs.)

- Quickie QM-720 BD:
No power Option Base w/ ASAP II Rehab Seat
- Quickie QM-720 SD:
Single Power Option Tilt (SPOT) or Recline base with ASAP II or Recline Rehab Seat
- Quickie QM-720 MPD:
Multiple Power Option Base (SPOT Tilt and Recline) with Recline Rehab Seat

Electronics Options:

VR-2 Non-Expandable Controls

- Industry standard non-expandable electronics
- Dependable, highly flexible, and simple to use
- 4-way toggle switch for operation of tilt and power legrests without the need for expandable controls

R-net Expandable Controls

- Assignable buttons for convenient control of common functions
- Optional full color display
- 8 uniquely programmable profiles for easy transitions between environments



QUICKIE® QM-7 Series

ASAP II Seat Frame:

- In-field width, depth, armrest, legrest, and back height adjustability
- Wide range of seat sizes to accommodate pediatric users
- Manual reclining backrest for increased function and sitting tolerance
- Optional fold down backrest for easy transport
- Extensive seating and positioning options
 - o Multiple armrest options
 - o Multiple legrest options
 - o Long thigh supports
 - o Fixed and articulating vent trays

Power Options:

- 50° of tilt to optimize positioning
- Power Recline with Power Shear Reduction
- 9" Powered Seat Elevator
- Power Elevating Legrests with Articulation
 - o Centermount
 - o Swing-Away



Quickie QM-7 Series Specifications

SPECIFICATIONS	QM-710 BC	QM-710 SC	QM-710 MPC	QM-715 HD BC	QM-715 HD SC	QM-720 BD	QM-720 SD	QM-720 MPD
HCPCS Code	K0848	K0856	K0861	K0850	K0858	K0868	K0877	K0884
Drive Wheel Position	Mid	Mid	Mid	Mid	Mid	Mid	Mid	Mid
Speed	up to 6.0 mph	up to 6.0 mph	up to 6.0 mph	5.0 mph	5.0 mph	8.0 mph	8.0 mph	8.0 mph
Overall Weight	310 lbs.	350 lbs.	375 lbs.	310 lbs.	350 lbs.	310 lbs.	350 lbs.	375 lbs.
Overall Width	25"	25"	25"	25"	25"	25"	25"	25"
Overall Length	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests
User Weight Capacity	300 lbs.	300 lbs.	300 lbs.	400 lbs.	400 lbs.	300 lbs.	300 lbs.	300 lbs.
Seat Widths	12" to 22"	12" to 22"	16" to 22"	16" to 22"	16" to 22"	12" to 22"	12" to 22"	16" to 22"
Seat Depths	12" to 22"	12" to 22"	16" to 20"	16" to 22"	16" to 22"	12" to 22"	12" to 22"	16" to 20"
Front Seat To Floor Height	17.5" to 20.5"	18" to 21"	16.5" to 19.5"	17.5" to 20.5"	18" to 21"	17.5" to 20.5"	18" to 21"	16.5" to 19.5"
Seat Frame Options	ASAP II	ASAP II or Recline Seat	Recline Seat	ASAP II	ASAP II	ASAP II	ASAP II or Recline Seat	Recline Seat
Power Seating Options	N/A	Tilt or Recline	Tilt & Recline	No Power Options	Tilt	No Power Options	Tilt or Recline	Tilt & Recline

COMPLEMENT YOUR QUICKIE



QUICKIE® Tension Adjustable Upholstery

Effective postural support and comfort in a lightweight, low profile design

E2611

JAY® Backrests



JAY® J3 Backrest

Flexible positioning and versatile hardware in an all-in-one solution


E2613, E2614, E2615, E2620



JAY® SPO Backrest

Exceptional positioning and postural management in a customizable package

E2613, E2617

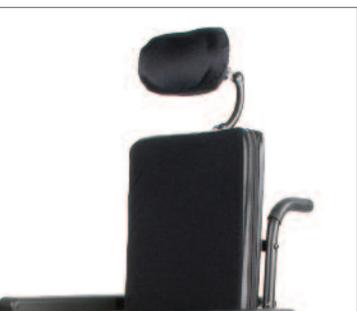

SPOT Curved Backrest

Solid support and industry-proven positioning components in a no-charge option

JAY® Cushions



JAY® J3 Cushion

Outstanding pressure distribution and optimal stability in a cushion that accommodates changing user needs

E2622, E2623, E2624, E2625



JAY® J2 Cushion

Exceptional skin protection and customization in a time-trusted cushion

E2622, E2623, E2624, E2625



JAY® Ion Cushion

Comfort, convenience and value in an accommodating skin protection and positioning cushion

E2607, E2608



JAY® Easy Cushion

Stability, skin protection and value in a foam and fluid cushion

E2607, E2608